


# Taleo<sup>®</sup> 1C50

Ready for everyday life.

Taleo was designed for active individuals who navigate varied indoor and outdoor environments and place a high value on effortless walking and the ability to go wherever life takes them.

**Designed for optimal outcomes:**

- Dual springs and a long carbon base enable a smooth rollover and efficient energy return at varying speeds
- The unique flexible connection of the carbon springs in the forefoot result in optimal adaption to varying ground conditions
- Customizable shock absorption with three different heel wedge options
- Protected against fresh, salt, and chlorine water and features water runoff channels on connection adapter and drainage holes in foot shell
- Cosmesis compatible with slim connection adapter

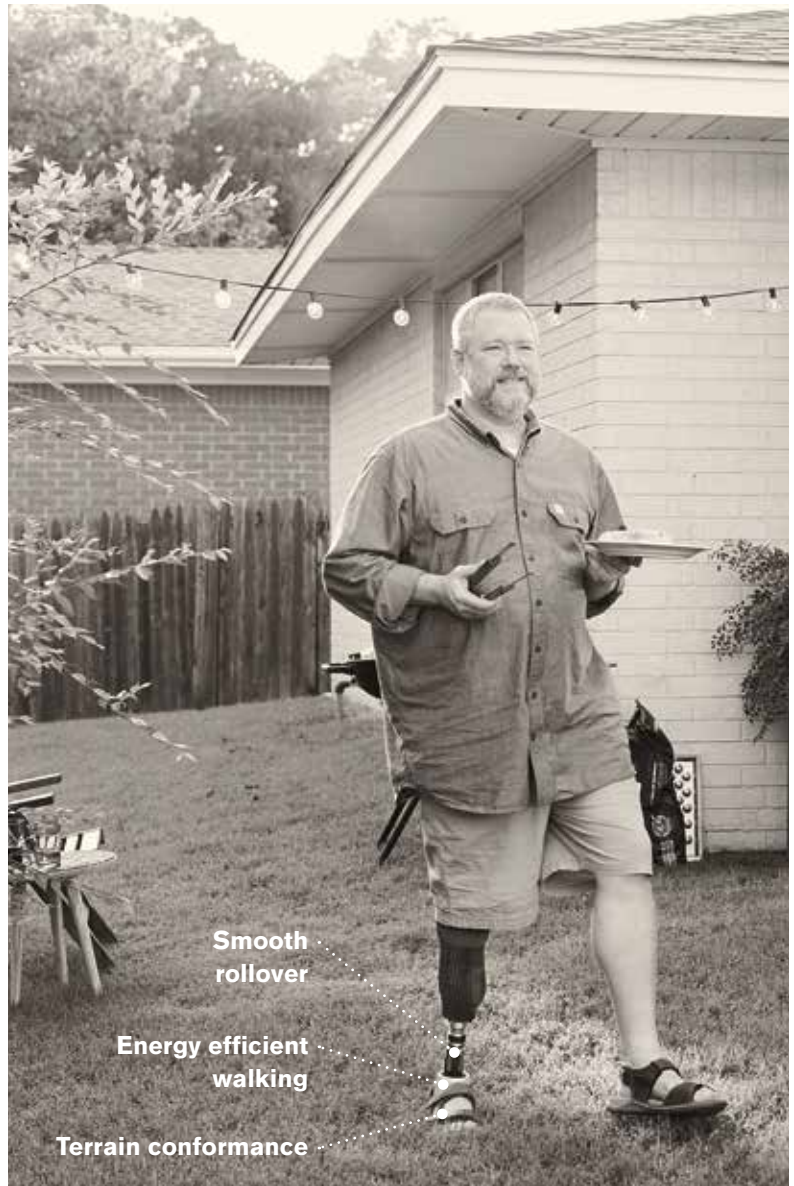

**Technical Data**

K Level	K3
Max. body weight	150 kg / 330 lbs
Sizes	22-30 cm
Footshell	Slim shape with 15 ± 5 mm heel height (size 22-25 cm) Normal shape with 10 ± 5 mm heel height (size 24-30 cm)
Weight without footshell*	450 g
Weight with normal footshell*	690 g
Structural height	150 mm (with normal footshell*)

\*Reference size is 26 cm.



# 1C50 Taleo



Order form for your Ottobock representative


Customer		Shipping address (if different from customer address)	
Customer no.	<input type="text"/>	Customer no.	<input type="text"/>
Company	<input type="text"/>	Company	<input type="text"/>
Street	<input type="text"/>	Street	<input type="text"/>
Postal code/city	<input type="text"/>	Postal code/city	<input type="text"/>
O&P professional	<input type="text"/>	Customer ref.	<input type="text"/>

## Selection of the spring stiffness relative to

1 body weight and activity as well as

Body weight lbs [kg]	Normal activity level	High activity level
Up to 113 [51]	1	2
114 – 129 [52 – 58]	2	3
130 – 149 [59 – 67]	3	4
150 – 170 [68 – 77]	4	5
171 – 195 [78 – 88]	5	6
196 – 220 [89 – 100]	6	7
221 – 254 [101 – 115]	7	8
255 – 287 [116 – 130]	8	9
288 – 330 [131 – 150]	9	–

2 foot size

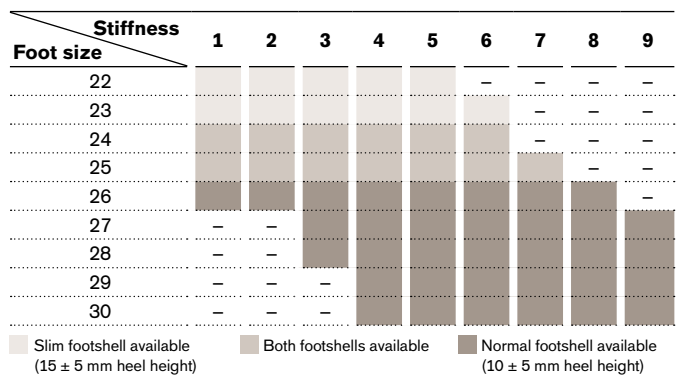







### Example for ordering 1C50=R26-4-P/4N

Article no.	Side	Size	- Stiffness	- P / Color	Shape
1C50=					
	Side	Size [cm]	Colour	Shape	
	Right R	22, 23, ..., 30	Beige 4	22–25 cm	S
	Left L		Light brown 15	Normal 24–30 cm	N



The scope of delivery includes the 1C50 foot module, 2C15 footshell including connection cover, 2F50 heel wedge set (comprising three different degrees of hardness) and a black Spectra sock.

### Recommended combinations\*

TT	TF
 <p>4R220 DVS + 6Y94= -F partial covered DVS liner and 453A30 ProFlex Plus knee sleeve</p>	 <p>7T451=1 SiOCX TF Pro; 6Y110 Skeo Sealing liner and 21Y21 ClickValve</p>
 <p>4R10=111 Quickchange</p>	 <p>3C98 C-Leg 3B1 Genium 3B5 Genium X3</p>
	 <p>4R57=WR Waterproof rotation adapter</p>
	 <p>3F1=1 + 99B120= Functional Cosmesis for C-Leg 4 3F1=2 + 99B120= Functional Cosmesis for Genium</p>

\* All components are sold separately and are available Ottobock products that are compatible with the 1C50 Taleo foot, which helps ensure optimal performance. Practitioners need to select components based upon individual patient criteria.